User manual

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Mode d'emploi

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Manual do utilizador

English	Ytango user manual	3
Nederlands	Ytango gebruikershandleiding	37
Deutsch	Ytango Gebrauchsanweisung	69
Français	Mode d'emploi du Ytango	101
Español	Ytango manual de usuario	133
Português	Ytango Manual do utilizador	165

Ytango user manual

Ytango

Congratulations on your purchase of the Ytango. The Ytango is the first behind-the-ear hearing aid that places both the microphone and speaker inside the ear canal for natural sound, directivity and reduced wind noise.

Table of contents

What is Ytango?	6	FDA Warning	25
Ytango's anatomy	9	Warnings	28
Ytango - module	10	Information about reparation services	30
Soft dome	12	Disclaimer	33
Battery	13	The data of your Ytango	34
The daily use of Ytango	14	Warranty	35
Place the battery	14	Manufacturer	197
Insert Ytango in the soft dome	16		
Place Ytango in the ear	17		
Remove Ytango from the ear	18		
Fitting procedure and configuration	18		
Switch between programs	19		
Telecoil functionality	21		
Maintenance and use	22		
Cleaning	22		
Storage	23		

What is Ytango?

The Ytango is a universal fit MaRiC™ (Microphone and Receiver in Canal) device using a soft dome.

The Ytango offers a range of important advantages:

- Natural sound. This is possible because the microphone is right in your ear, so you hear with natural directivity.
- Clear sound quality by positioning of the speaker close to the eardrum.
- Almost invisible.
- Open fit dome possible for minimal occlusion (plugged ear feeling, and/or distortion of your own voice).
- Minimum wind noise due to the microphone positioned deep in the ear
- Soft dome design protects the instrument against damage caused by earwax.
- Suitable for cell phone use.

- Suitable for sporting activities.
- Easy to use and maintain.
- Contemporary design.

The Ytango product includes the following features:

- An adaptive feedback cancellation system - this prevents spontaneous 'whistling' of the device which could happen as you insert it in your ear. This is normal during the insertion process. As soon as you have placed it in your ear canal, the whistling will stop and the adaptive feedback system will protect you from feedback.
- An adaptive noise reduction system - this improves the hearing comfort in noisy environments.
- Data logging; a feature that stores daily device interactions and can be analyzed during examinations.

In addition the Ytango Pro and Ytango Pro-T add the following features:

- Four different memories - Use up to four different listening situation settings when programming the Ytango. You may request that your dispenser programs the memories to your special environmental needs. When the Ytango has been set up with more than one memory, the memories can be selected by AirTAP functionality. For a description about this functionality, please refer to the chapter 'Switching between programs' (page 19).
- The Ytango Pro-T features a telecoil connectivity (T-coil) which allows audio sources to be directly connected to the hearing aid from devices such as special telephones, FM systems (with neck loops), and induction loop systems (also called "hearing loops").

	Ytango	Ytango Pro	Ytango Pro-T
MaRiC® for natural hearing	•	•	•
Natural directionality	•	•	•
AirTAP™ (4 programs)		•	•
Telecoil			•
Soft dome selection	Open (various sizes) Closed (various sizes)		
Colors	Black, Titanium, Champagne		
Fitting range	60 dB (open) 75 dB (closed)		

Ytango's anatomy

A Ytango - module


- 1 Battery compartment
- 2 Pre-bent wire
- 3 MaRiC unit
- 4 Microphone opening
- 5 Speaker opening
- 6 Telecoil switch

B Soft domes

- 7 Sound opening
- 8 Open dome
- 9 Closed dome

C Battery (type 312)

- 10 Plus pole (+)
- 11 Minus pole (-)

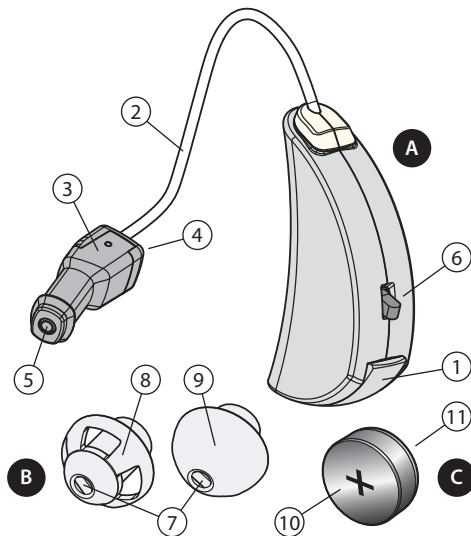Ytango's anatomy

Ytango - module

The module uses an advanced programmable microprocessor to ensure the processing and amplification of the sound is adapted to your particular hearing loss. The casing of Ytango module has been made of durable plastic. For the maintenance of the casing, please refer to the chapter 'Maintenance and use' (page 22).

Battery compartment

This compartment houses the battery that provides the power for your Ytango. The Ytango can be switched off by slightly opening the battery door to the first click. Open the door completely to replace the battery.

Right and Left MaRiC cables

Right and left MaRiC cables can be identified with their color coding:

- Red: right ear
- Blue: left ear

The MaRiC connector cables are available in various sizes. It is important that the right size is selected in order to avoid that the MaRiC unit is positioned incorrectly in the ear. A correctly positioned cable positions the MaRiC unit in the middle of the ear canal. Your dispenser will help determine which MaRiC cable size is suitable for your ear.

Ytango's anatomy

MaRiC unit (Microphone and Receiver in Canal)


The Ytango uses your own outer ear to capture and focus sound, just as nature intended. The microphone is placed inside the ear, to pick up and amplify sound waves with natural and clear quality.

Microphone opening

The microphone picks up the sound, which is processed and amplified by the Ytango.

Speaker (or receiver) opening

The processed and amplified sound emits through speaker opening and is directed via the soft dome towards the ear drum.

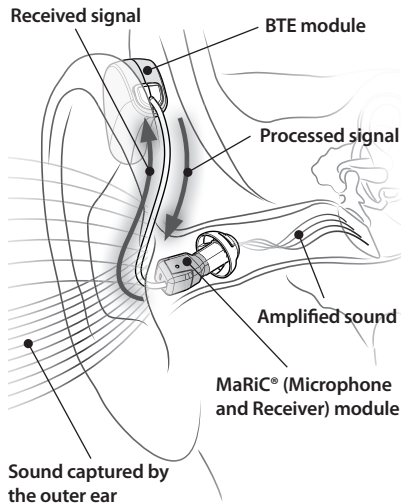Ytango's anatomy

Soft dome

The Ytango uses soft domes available in various sizes and shapes. The size and shape of the soft dome to be used depends on the size and shape of your ear canal. Choosing the correct size soft dome for your Ytango is extremely important for comfort and functionality.

Oversized soft dome may cause discomfort and occlusion. Undersized soft dome may lead to MaRiC movement in the ear canal causing feedback (or whistling). Your dispenser will determine the soft dome that suits you best during your fitting.

The soft dome of the Ytango is made of silicone and can be replaced when necessary. One of the most common hearing instrument repairs is for a damaged or clogged receiver due to earwax. The Ytango soft domes have a

built-in wax protection system that prevents the intrusion of earwax into the device. It is therefore important to use only the soft domes that have been designed specifically for the Ytango. For its maintenance and cleaning, please refer to the chapter 'Maintenance and use' (page 22). The illustration may be different from your soft dome.

Sound opening

The processed and amplified sound emits through this opening and is directed towards the eardrum.

Ytango's anatomy

Open fit soft domes

Several openings (holes) can be found on the dome. Generally speaking, open soft domes will have less amplification in the lower frequencies than a more closed soft-dome. Open-fit soft domes are better suited for a mild to moderate hearing loss.

Closed fit soft domes

Closed fit soft domes support a wider range of hearing losses. A side effect of closed-fit soft domes is that your own voice will sound louder than normal - this is called occlusion.

Battery


Ytango uses type 312 batteries indicated with a brown sticker. The device is switched off upon slightly opening the battery door to the first click. Replacing a battery is possible with the battery compartment door completely opened.


The daily use of Ytango

Before you start any of the actions mentioned in chapters 'The daily use of Ytango' and 'Maintenance and use', make sure you are using a clean, horizontal working surface.

Inserting the battery

Ytango uses type 312 batteries indicated with a brown sticker. Average battery life is around 150 to 180 hours of continuous use. It is advised to switch the hearing aid off when it is not used, to prolong battery life. The device is switched off upon slightly opening the battery door to the first click. Replacing a battery is possible with the battery compartment door completely opened.

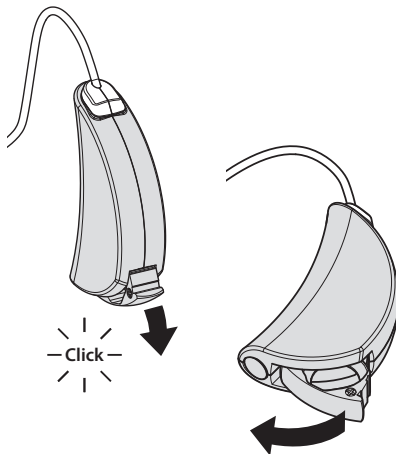
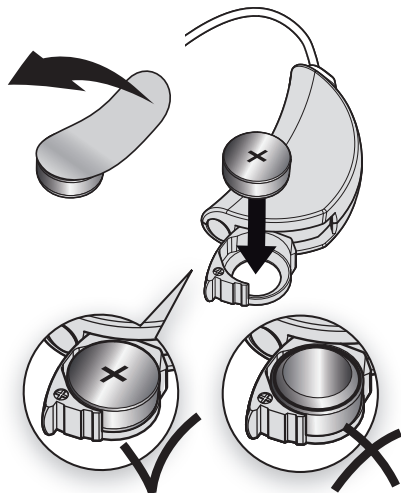
To insert a new battery: first remove the brown sticker from the battery's "+" pole. Open the Ytango battery door and insert the battery as indicated by the battery door design and "+" marking.

The daily use of Ytango

Insert Ytango's MaRiC in the soft dome


Take off the old dome using your fingernails. Push the new dome onto the speaker. Give it a push to secure it safely.

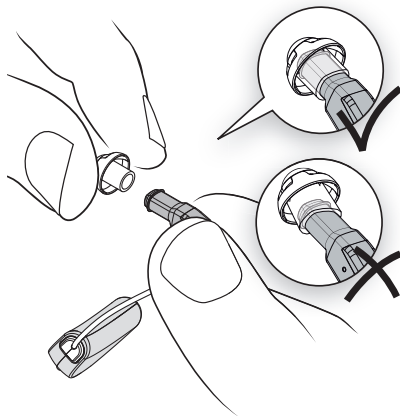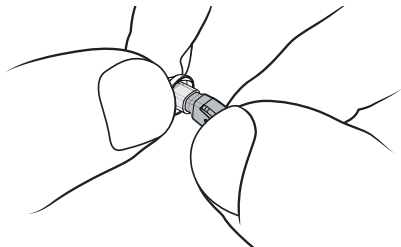Placing the Ytango MaRiC in the ear

Left and right devices can be identified with their color coding on the MaRiC cable:

- Red: right ear
- Blue: left ear
- Place the MaRiC with soft dome loosely in the ear.
- Insert the device into the ear by pressing on the outer end of the MaRiC. Continue to gently push the device deeper into the ear. If your mouth is opened wide, insertion will be easier.

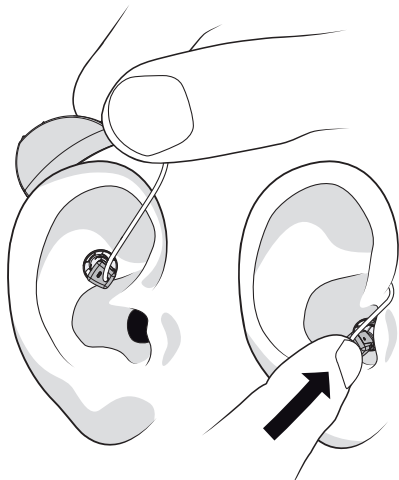The daily use of Ytango

Inserting a new soft dome into the ear can be simplified by using some ear lubricant or crème (a tiny amount is sufficient).

Attention: Be careful to make sure the ear lubricant does not clog the sound opening of the soft dome!

Ideally, the flat side of the module (where the microphone is situated) will be flush to the entry of the ear canal after insertion. In this way, the device is practically invisible and the sound optimal.

Make sure your hands are clean while inserting the Ytango. Make-up, hand lotion, body oils, hair gel, hairspray and similar products could make the module dirty or contribute to clogging of the microphone opening.

Removing Ytango from the ear

Hold the cable firmly and pull it slowly out of the ear. Remove Ytango module beyond your ear.

Fitting procedure and configuration

The fitting procedure and configuration of the hearing aid is both an important and individual process. A specialized computer program with specific connectors are available to your dispenser to accomplish this. For this reason, consult your dispenser for any fitting or configuration changes. Your dispenser is the authorized person to carry out any necessary modifications.

Switch between programs (Ytango Pro and Ytango Pro-T only)

If you are wearing a Ytango Pro or Ytango Pro-T and your instrument has been configured with more than one memory, memories can be selected using AirTAP™ functionality which allows you to choose memory for a different environment.


To change between the memories, please gently tap your ear with hand. Most people require a bit of practice to figure out the right method. General guidelines are: tap with a closed “cupped” hand (fingers together) gently. A tone will be heard indicating the number of the new memory that was entered.

The daily use of Ytango

The Ytango Pro AirTAP is defaulted to use a single tap to change memories. This mode requires the user to tap the ear once to switch memories. The dispenser can change this to a two-tap memory change if desired. This mode requires the user to tap the ear twice to switch to the next memory. Should the user not want to use the different programs, the AirTAP function can also be disabled.

After inserting the battery, the Ytango instrument will start in program one. Gently tapping the ear will produce two beeps, indicating that program two has been activated. After another gentle tap, three beeps will be heard indicating that program three has been activated. A maximum of four programs can be set up. While program four is active, a gentle tap sets the device back into memory one. One beep will be heard indicating memory one is active.

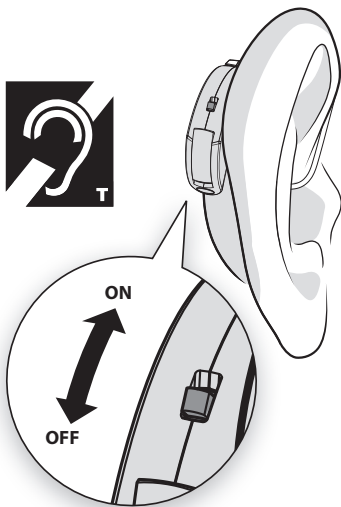Telecoil functionality (Ytango Pro-T only)

The telecoil functionality is controlled by the switch located on the back of the Ytango. The telecoil is inactive when the switch is in down position. To activate the Telecoil place the switch in the up position. The acoustic microphone will be shut off when using the telecoil feature.

Maintenance and use

The build up of earwax differs from person to person. In case you observe that the sound outlet of the soft dome has been clogged, you may clean the soft dome (see paragraph 'Cleaning') or replace it.

Check the module(s) and the soft dome(s) daily for clogs and damages. This guarantees you optimal performance. If you notice that the soft dome is damaged, it is important that you do not continue to use it. Replace it immediately! To keep the performance of your device optimal and minimize the chance of infection, we advise you periodically replace the soft dome with a new one. Soft domes are available at your dispenser.

Cleaning

MaRiC (Microphone and Receiver in Canal) module:

- Clean daily with a soft and dry tissue.
- In case of severe dirt, wipe with special (moist) cleaning tissues for hearing aids. Be careful to avoid moisture entering the small holes for the microphone or the speaker . This could damage the MaRiC

Ytango (behind the ear) module:

- Clean daily with a soft and dry tissue.
- In case of severe dirt, you may clean it with special (moist) cleaning tissues.

Soft dome (or soft tip):

- Attention: soft domes should be cleaned separately from the MaRiC module. The MaRiC module could be damaged if the softdomes are cleaned together with the MaRiC.
- Clean daily with a soft and dry tissue.
- In case of severe dirt, you may clean it with special (moist) cleaning tissues or with soluble cleaning tablets for hearing aid products. You may alternatively have it cleaned with an ultrasound cleaner.
- Clogged soft domes can be rinsed off under streaming lukewarm water. In case of severe dirt, change in the properties or damage of the material occurs, the soft dome should be replaced.

Attention: Never use aggressive cleaning products, abrasives or solvents! This can cause damage to the soft dome or the module.

Storage

At night:

Open the battery door slightly to shut the device off. This will increase the operational life of the battery. Store the device in a secure location outside of its box. This allows for the evaporation of any moisture in the device.

Maintenance and use

Short-term storage:

Use the provided carrying box when you want to remove the device from the ear temporarily.

Attention: do not carry the hearing aid unprotected in your pocket, purse, or bag.

Long-term storage:

When not using your hearing aid for longer periods of time, clean the module and soft dome before storing them inside the carrying box. Always remove the battery from the device in this situation.

FDA Warning

The following additional information is provided in compliance with U.S. Food and Drug Administration (FDA) regulations:

Warning to hearing aid dispensers

A hearing aid dispenser should advise a prospective hearing aid user to consult promptly with a licensed physician (preferably an ear specialist) before dispensing a hearing aid if the hearing aid dispenser determines through inquiry, actual observation, or review of any other available information concerning the prospective user, that the prospective user has any of the following conditions:

i. Visible congenital or traumatic deformity of the ear.

ii. History of active drainage from the ear within the previous 90 days.

iii. History of sudden or rapidly progressive hearing loss within the previous 90 days.

iv. Acute or chronic dizziness.

v. Unilateral hearing loss of sudden or recent onset within the previous 90 days.

vi. Audiometric air-bone gap equal to or greater than 15 decibels at 500 Hertz (Hz), 1,000 Hz and 2,000 Hz.

vii. Visible evidence of significant cerumen accumulation or a foreign body in the ear canal.

viii. Pain or discomfort in the ear.

RF Immunity Level

These hearing instruments have a cell phone immunity rating of M4. For your cell phone to be compatible with these hearing instruments, the cell phone needs an immunity rating of

FDA Warning

M1 or higher. Please consult your cell phone specifications for your cell phone immunity rating.

Important notice for prospective hearing aid users

Good health practice requires that a person with a hearing loss have a medical evaluation by a licensed physician (preferably a physician who specializes in diseases of the ear) before purchasing a hearing aid. Licensed physicians who specialize in diseases of the ear are often referred to as otolaryngologists, otologists, or otorhynolaryngologists. The purpose of the medical evaluation is to assure that all medically treatable conditions that may affect

hearing are identified and treated before the hearing aid is purchased.

Following the medical evaluation, the physician will give you a written statement that states that your hearing loss has been medically evaluated and that you may be considered a candidate for a hearing aid. The physician will refer you to an audiologist or hearing aid dispenser, as appropriate, for a hearing aid evaluation. The audiologist or hearing aid dispenser will conduct a hearing aid evaluation to assess your ability to hear with and without a hearing aid. The hearing aid evaluation will enable the audiologist or dispenser to select and fit a hearing aid to your individual needs.

If you have reservations about your ability to adapt to amplification, you should inquire about the availability of a trial-rental or

purchase-option program. Many hearing aid dispensers now offer programs that permit you to wear a hearing aid for a period of time for a nominal fee after which you may decide if you want to purchase the hearing aid.

Federal law restricts the sale of hearing aids to those individuals who have obtained a medical evaluation from a licensed physician. Federal law permits a fully informed adult to sign a waiver statement declining the medical evaluation for religious or personal beliefs that preclude consultation with a physician. The exercise of such a waiver is not in your best health interest and its use is strongly discouraged.

A hearing aid will not restore normal hearing and will not prevent or improve a hearing

impairment resulting from organic conditions. In most cases infrequent use of a hearing aid does not permit a user to attain full benefit from it.

Children with hearing loss

In addition to seeing a physician for a medical evaluation, a child with a hearing loss should be directed to an audiologist for evaluation and rehabilitation since hearing loss may cause problems in language development and the educational and social growth of a child. An audiologist is qualified by training and experience to assist in the evaluation and rehabilitation of a child with a hearing loss.

Warnings

The soft dome is an essential part of the Ytango. Never use the product without a soft dome.

Once the device has been fitted to one specific ear, you should not wear it in the other ear. This may lead to a deterioration of hearing performance or hearing loss.

Never use damaged or dirty soft domes. This can negatively affect the hearing performance or result in an uncomfortable fit.

Replace the soft dome of your device regularly, even if it is not visibly damaged or dirtied, see chapter 'Maintenance and use' (page 22). When using Ytango for the first time, begin with a continuous wearing period of 4 hours and subsequently increase wearing time with one hour a day. If you suffer from pain,

continuous irritation or itching, do not use the devices before you receive specialist advice. If your complaints persist, consult a physician. Have your ears and your hearing checked regularly, both during the beginning period and afterwards.

If for some reason it is not possible to remove Ytango from your ear, consult your dispenser or a physician.

Do not use the Ytango in the following cases:

- A medical history of the ear.
- Exostoses (bulging bone growth) in the ear canal.
- Heavy earwax secretion.

The device should not be:

- dropped.
- exposed to moisture or water.

- immersed in water.
- used while swimming, taking shower, bathing or in the sauna.
- cleaned wet.
- exposed to extreme high temperatures, such as when using a hair dryer.
- worn in mines or other environments where explosion risk is high.
- exposed to X-ray radiation or strong magnetic fields, such as during medical investigation.
- worn in the places where hairspray (hairdresser) or solvents (such as dye) are used.

Care should be taken in using make-up, crème, body oil, gel, hairspray and similar products, while you are wearing or inserting the Ytango. They can dirty up your device and, for example, contribute to the clogging of the microphone opening.

Keep the Ytango away from small children and pets. In case of swallowing one of the components, there is a choking hazard to small children.

If a battery leaks and its content touches the skin, flush exposed skin with copious quantities of flowing lukewarm water for a minimum of 15 minutes. Wash skin with soap and water.

If a battery is swallowed, immediately go to the hospital or see a doctor.

Used batteries and old hearing aids should be collected and disposed as chemical waste.

Information about reparation services

The Ytango has been developed by ExSilent BV, Overschiestraat 63, Amsterdam. For your complaints or reparation requests, you should contact the provider from whom you purchased your device. Complaints and reparations settlement may vary by provider.

Limited Manufacturer's Warranty.

This limited guarantee will be applicable in any part of the world. ExSilent BV ("ExSilent") guarantees that this ExSilent product ("Product") is free from defects in materials, design and finish at the time of its original purchase subject to the following provisions and conditions:

1. This limited guarantee is given to the end user of the Product ("Client"). It is given without exclusion or restriction of (i) any of the Client's binding statutory rights or (ii) any of the Client's rights against the seller/concessionaire of the Product.

2. The guarantee period is twenty four (24) months from the date on which the first Client purchased the Product. In the event of a successive purchase or other change of owner/user, the guarantee period will continue for the remaining part of the twenty four (24)-month period and remain unaffected in all other respects.

3. During the guarantee period, ExSilent or its approved service provider will repair or replace the defective Product (fully at the discretion of ExSilent). ExSilent will return the repaired Product or another Product to the Client in good condition. All components or other equipment which are replaced will become the property of ExSilent.

4. A replaced product will be given an extended or renewed guarantee period. In case of a repaired product, the old guarantee will continue to apply.

5. This limited guarantee does not cover painted covers or other comparable personalized components.

6. This limited guarantee is not applicable to damage caused by normal wear and use. Neither will this limited guarantee be applicable if: (i) the defect was caused by the fact that the Product was subjected to: use in conflict with the owner's or user directions, rough treatment, exposure to moisture, damp or extreme thermal or environmental conditions or rapid changes in such conditions, corrosion, oxidation, non-permissible alterations or connections, non-permissible openings or repairs, repairs made with non-permissible replacement parts, incorrect use, improper installation, accidents, force majeure, spilling of food or liquids, the influence of chemical products or other facts over which ExSilent cannot reasonably have

control (including but not limited to defects in consumable components, such as batteries which only have a limited life by nature) unless the defect has been caused directly by defects in materials, design or finish; (ii) ExSilent or its approved service provider received no notification of the defect from the Client within thirty (30) days of the occurrence of the defect within the guarantee period; (iii) the Product was not returned to ExSilent or its approved service provider within thirty (30) days of the occurrence of the defect within the guarantee period; (iv) the serial number was removed, erased, damaged, changed or is illegible; (v) the defect was caused by the fact that the Product was used with or was connected to an additional component which was not manufactured or supplied by ExSilent or used for a purpose other than that for which it was intended; (vi) The defect was caused by the fact

Information about reparation services

that the battery was subjected to a short-circuit or that the connecting pieces of the battery housing or cells were broken or show evidence of tampering or that the battery was used in another device than that for which it is intended.

7. In order to make a claim under the limited guarantee, the Client must send a legible and unaltered, fully-completed original guarantee card to ExSilent or its authorized service provider within two weeks of purchasing the Product.

8. This limited guarantee is the sole and exclusive remedy available to the Client against ExSilent and the sole and exclusive liability of ExSilent towards the Client for defects in or poor functioning of the Product. This limited guarantee replaces all other guarantees and liabilities, including verbal, written, (non binding) statutory, contractual, unlawful and

otherwise. ExSilent is not liable for any indirect, secondary or consequential damage, costs or expenses. Neither is ExSilent liable for any direct damage, costs or expenses, if the Client is a legal entity.

9. Any amendment or change to this limited guarantee requires the prior written approval of ExSilent.

Disclaimer

© 2011 ExSilent. All rights reserved.

Illegal reproduction, transfer, distribution or holding of this document, or a part hereof, in one way or another, without the prior written consent of ExSilent is forbidden.

ExSilent reserves the right to change and make improvements in its products that are described in this document without prior notification.

In no case, ExSilent can be held liable for any loss of data, neither for income, nor for any special incidental, nor for indirect damage.

The content of this document is provided without any warranty. Unless the applicable law requires otherwise, neither explicitly nor implicitly nor guarantee will be given in regard to the exactitude, reliability and content of this document. This includes but is not limited to implicit warranties regarding its marketability and suitability to attain a certain goal.

ExSilent reserves the right at any time to change and to revoke this document without prior notification.


The data of your Ytango

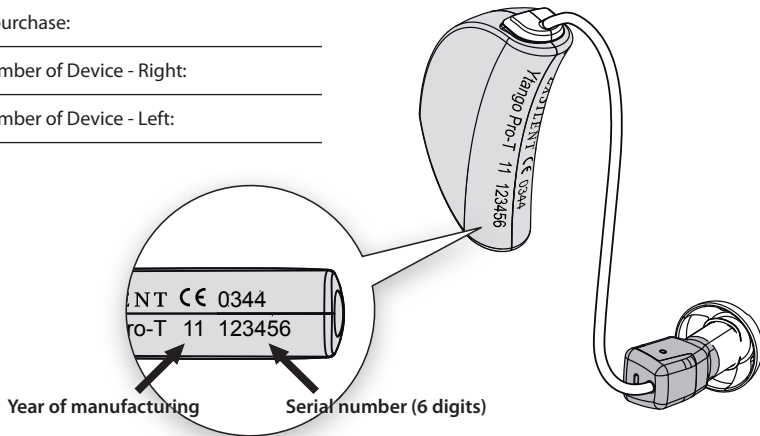Please note the data of your Ytango below.

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ExSilent B.V. (legal manufacturer)